


Digital baby scale with fine graduation

seca 717 




▲ Fine 2 g graduation to determine the smallest changes in weight . HOLD and TARE functions for ease of operation . Automatic switch-off reduces power consumption . Robust and easy to clean weighing tray . Can be fitted with optional measuring rod seca 231

A safe and secure start to life

A new baby has arrived – parents and doctors are there to protect it. Every gram is important. With its ease of operation and fine 2 g graduation, the digital baby scale seca 717 is a valuable aid when controlling the development of a newborn child, especially in hospital. For healthy growth from the very beginning.

Fine graduation for determining the smallest changes in weight.

The sensitive electronics of the seca 717 determine the weight of the newborn child down to 2 grams exactly. As far as precision is concerned, this scale is far ahead of its competitors. The intelligent damping system ensures that the measuring result is still precise even if the baby should become restless. The result is easy to read on the brightly illuminated LED display even if there is little light.

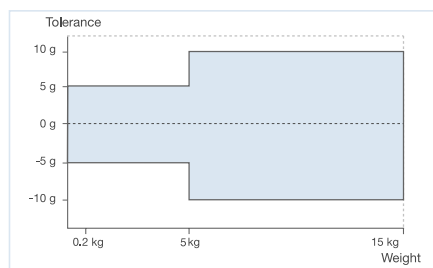

TARE and HOLD functions for ease of operation.

With the integrated TARE function, the weight of a soft mat or nappy can be easily deducted at the touch of a button. With the HOLD function, the determined weight can be stored until the baby has been attended to. As the zero point of the scale can change as the result of external influences especially in the case of such highly sensitive electronics, it is automatically checked and aligned every second. The automatic switch-off after ten minutes prevents unnecessary power consumption.

Robust and easy to clean weighing tray.

The baby can easily be placed on the absolutely tilt-resistant weighing tray of the seca 717, the sides of which are open. The baby feels safe and secure in the ergonomically shaped, slightly curved tray. The surface material is kind to the skin and is easy to clean and disinfect.



Optional measuring rod seca 231

When fitted with the optional measuring rod seca 231, weighing and measuring can be carried out practically and easily in one step. This saves the baby a lot of unrest and the nursing staff a lot of time and work.

Technical data

seca 717

Capacity: 15 kg

Graduation: 2 g < 6 kg > 5 g

Approval class: III

Dimensions (WxHxD): 551 x 163 x 320 mm

Weight: 6.5 kg

Power supply: Mains adapter

Functions: TARE, HOLD, weighing range switch-over

Optional: Measuring rod seca 231, head and foot positioner seca 429, transport case seca 425

